

# iBioSim Grid Tutorial

Jason Stevens

August 25th, 2011

## Contents

<b>1</b>	<b>Introduction</b>	<b>2</b>
<b>2</b>	<b>Creating a Grid Component</b>	<b>2</b>
<b>3</b>	<b>Creating and Populating a Grid</b>	<b>3</b>
<b>4</b>	<b>Simulating and Visualizing a Grid</b>	<b>5</b>


# 1 Introduction

This tutorial will step you through the process of creating a spatial grid model of a cell population using iBioSim. Before following this tutorial, please complete the basic iBioSim tutorial and read through the manual. Knowledge of some of the material discussed in those resources will be assumed.

## 2 Creating a Grid Component


A grid is a structure used for separating components spatially. So the first thing to do is create the components that will be placed onto the grid. This tutorial is going to use a multicellular circuit that will create a spatial NOR function. (For a similar NOR circuit that has been built and characterized, see Tamsir et. al 2011. For a collection of multicellular logic gates, see Regot et al. 2011.) A NOR gate gives an ON response when both inputs are OFF; otherwise the circuit is OFF. This circuit is comprised of three different cells types: two sender cells that create diffusible input species, and a cell that receives the inputs and makes the NOR computation.

1. Below is the Sender 1 model. It contains two constitutive species, one which will act as an input (SignalA) and the other as an aesthetic reporter (CFP). Make sure that SignalA is marked as diffusible and that the model is enclosed in a compartment (tick the box on the toolbar).

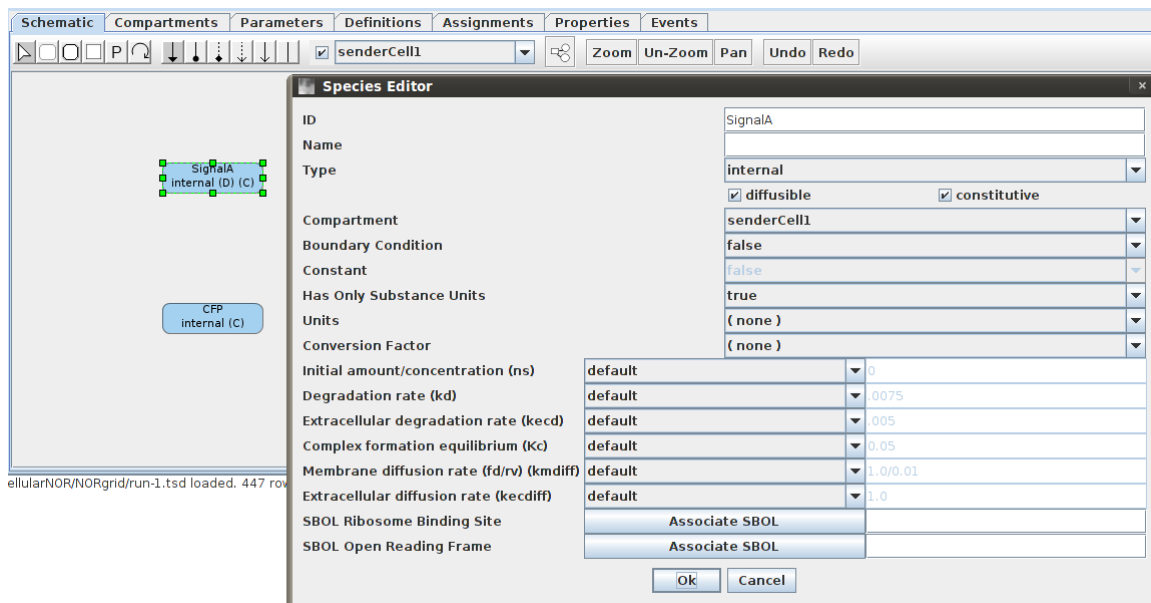

2. For the Sender 2 model, create a model similar to Sender 1 that has a different name (e.g., SignalB) for the diffusible species.
3. Below is the model for the NOR gate cell. At the left, you can see that two species repress a promoter that leads to production of a GFP species. If both of the repressing species are marked diffusible and given the same ID as the species in Sender 1/2 (e.g., SignalA/B), then membrane diffusion reactions will be created and any SignalA/B species can come into or diffuse out of the cell. If these species were not present or not marked diffusible, then no membrane diffusion reactions would be created. In the figure there is also a species panel


which shows custom diffusion rates that should be used for both SignalA and SignalB. Also, please change Kr to 5 for the two repression influences (by double-clicking on the red repression influences) and ko to 0.025 for the promoter. Like Sender 1/2, this model needs to be enclosed in a compartment. Enclosing it in a compartment creates a virtual membrane, which species can be transported through via membrane diffusion reactions. (Note that membrane diffusion rates are taken from the inner species. The membrane diffusion rates for the Sender models are kept as default for the sake of simplicity.)

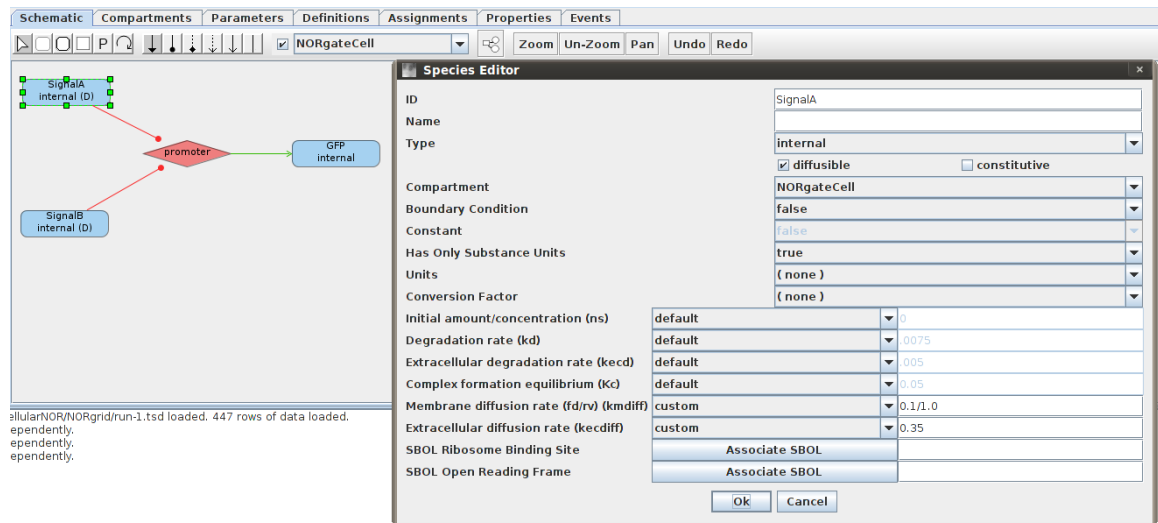


4. The grid components are now done, so let's put them onto a grid!

### 3 Creating and Populating a Grid


1. To create a grid, create a new model and tick the "Make Grid" box when entering the model ID. You will then be prompted to create a new grid. Make this grid 6 x 15 and fill it with the NOR gate model. After a moment, you will see a populated grid in the schematic. If the grid is too large to fit, you can zoom out by selecting the zoom button in the toolbar, or holding alt or control, and then using the mouse wheel (or whatever you use for vertical scroll).

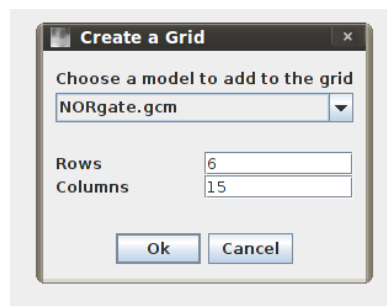


2. The grid now only contains the NOR gate model, so we need to add the Sender models to the grid as well. To do this, left-click on the four corners of the grid (they should appear light-blue if selected), then right-click and choose the "Clear Selected Location(s)" option. (Placing the Sender models at the corners will allow for a good demonstration of the NOR functionality.)

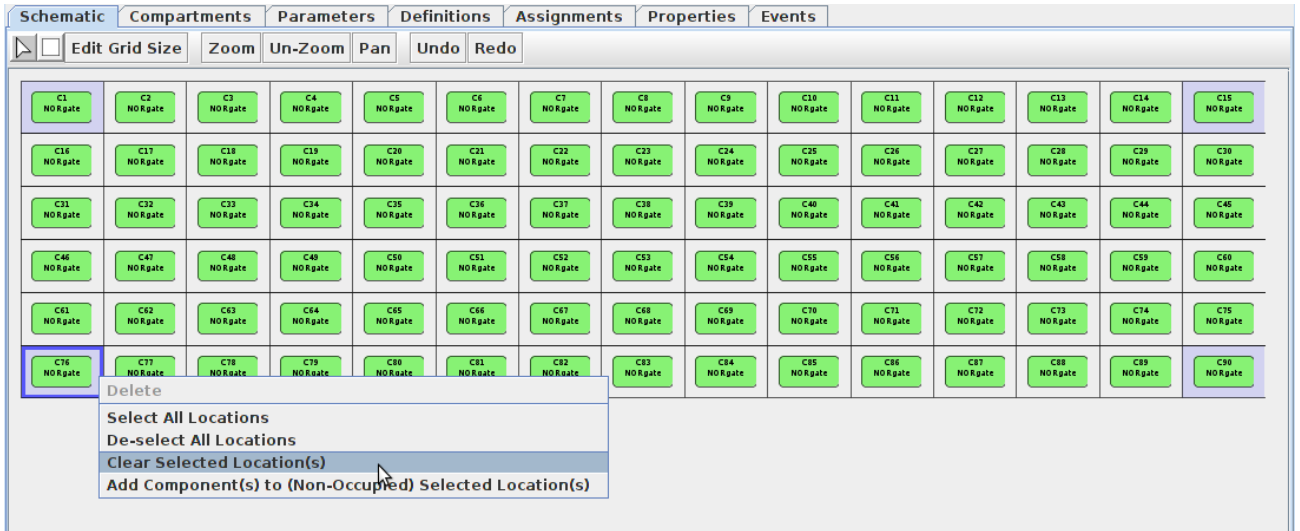

- Now let's add Sender 1 to the top two corners. De-select the bottom two corners (by single-clicking), then right-click and choose the "Add Component(s) to Non-Occupied Selected Location(s)" option. Choose the Sender 1 model at the prompt.

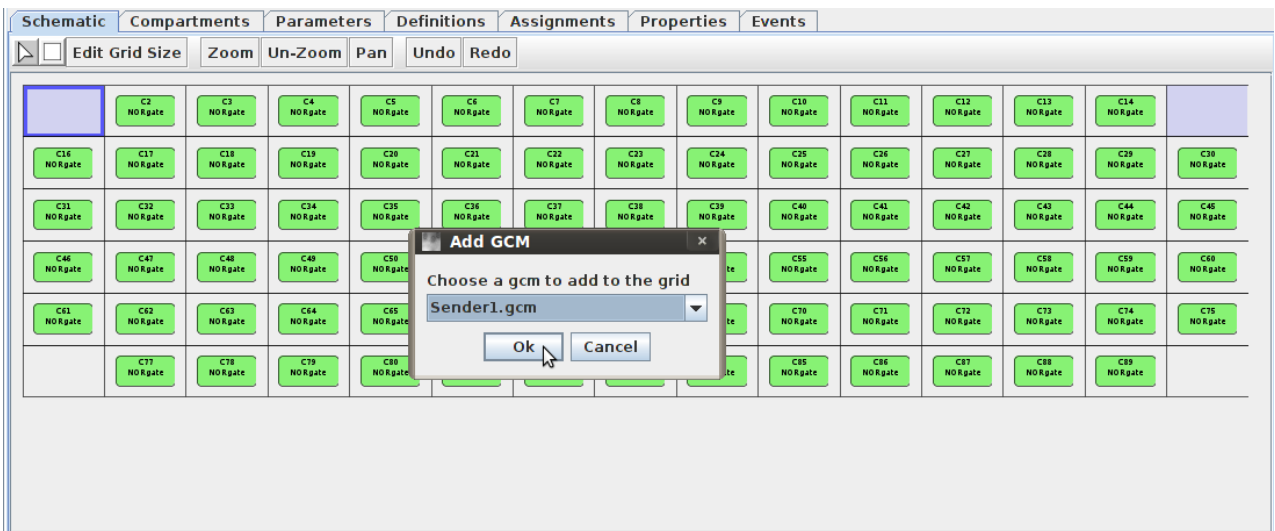


- Now let's add Sender 2 to the bottom two corners. To demonstrate another way of adding components, select the component radio button on the toolbar (looks like a square). Left-click on one of the empty bottom two corners then choose the Sender 2 model at the prompt. Do the same for the remaining empty corner.

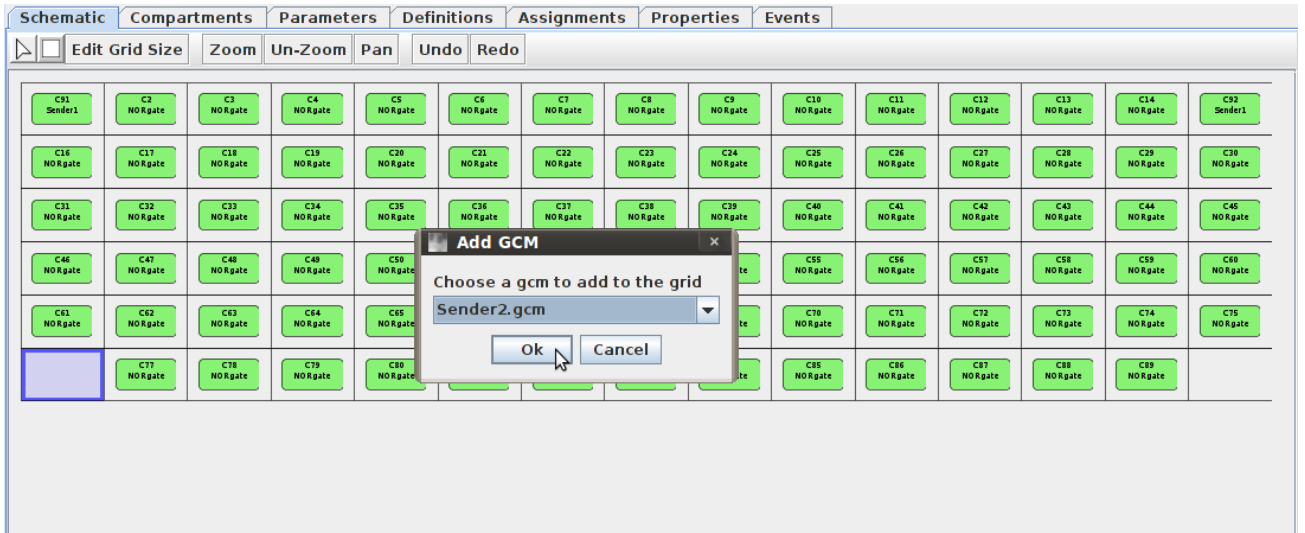

5. The grid is now complete, so let's simulate it!

## 4 Simulating and Visualizing a Grid

1. Simulating a grid model is no different than simulating a non-grid model. For my simulation run, I turned on abstraction, used Monte Carlo/Gillespie, and set a time limit of 400. When you've chosen your simulation options, click the little green play button on the top toolbar. This will generate the SBML and then simulate the model. Simulation should take under a minute with the options provided.

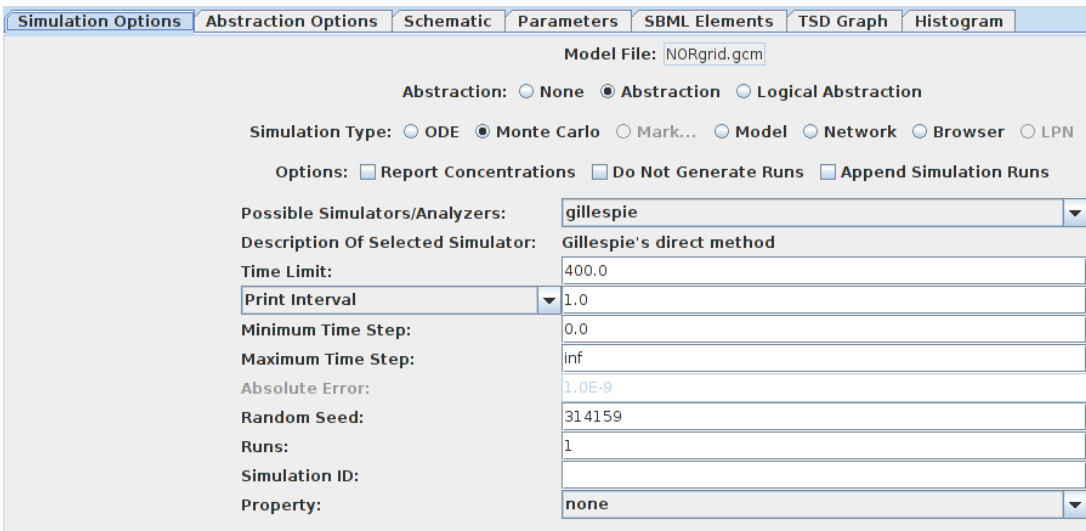

2. When the simulation has finished, click on the Schematic tab. If the model is too large to fit on your screen, you can zoom out by holding control or alt and using the mouse wheel (or whatever you use for vertical scroll). Now click the "Choose Simulation" button and choose the simulation run you would like to use for playback (probably "Run-1.tsd").
3. Now it's time to choose appearance schemes to maximize our ability to see what's happening during simulation playback. To alter a species' appearance, select its compartment or grid

The screenshot shows the 'Simulation Options' window with the 'Abstraction Options' tab selected. A dialog box titled 'Select inner species appearance scheme' is open, displaying options for the 'GFP' species. The dialog includes fields for 'Species' (GFP), 'Color Gradient' (green), 'Size Gradient' (none), 'Opacity Gradient' (none), 'Min. number of molecules visible' (0), 'Saturating number of molecules' (10), and 'Apply to' (all components with this model). 'Ok' and 'Cancel' buttons are at the bottom.

- 6



5. Hooray! You have now completed the grid tutorial. Hopefully this has given you a good sense of how to use iBioSim for static spatial modeling. If you have any questions, comments, suggestions, or issues, you can send an email to [atacs-bugs@vlsigroup.ece.utah.edu](mailto:atacs-bugs@vlsigroup.ece.utah.edu). Be sure to include "BUG", "CHANGE", or "FEATURE" (even if it's a question) in the subject line or it will be marked as spam. Happy modeling!