

iBioSim Installation Instructions

Chris J. Myers

August 26, 2011

Contents

1	General Requirements	2
2	Installation on Windows	2
3	Installation on Linux	4
4	Installation on MacOS	8

1 General Requirements

There are versions of iBioSim available for Windows, Linux, and MacOS. You can download the appropriate installation file from:

<http://www.async.ece.utah.edu/iBioSim>

iBioSim requires that you have Java Runtime Environment 1.5 or higher install on your system. It is also useful to have Graphviz. Also, you should associate files with a “.dot” extension with the Graphviz tool, and files with the “.xhtml” extension to your web browser.


2 Installation on Windows

Download and execute iBioSim-**<version>**-Setup.exe.


The installation uses InstallJammer. It first asks you for your preferred installation language. Make your selection and press OK.




The next screen tells you that you what version you are installing. Press **Next** to continue.

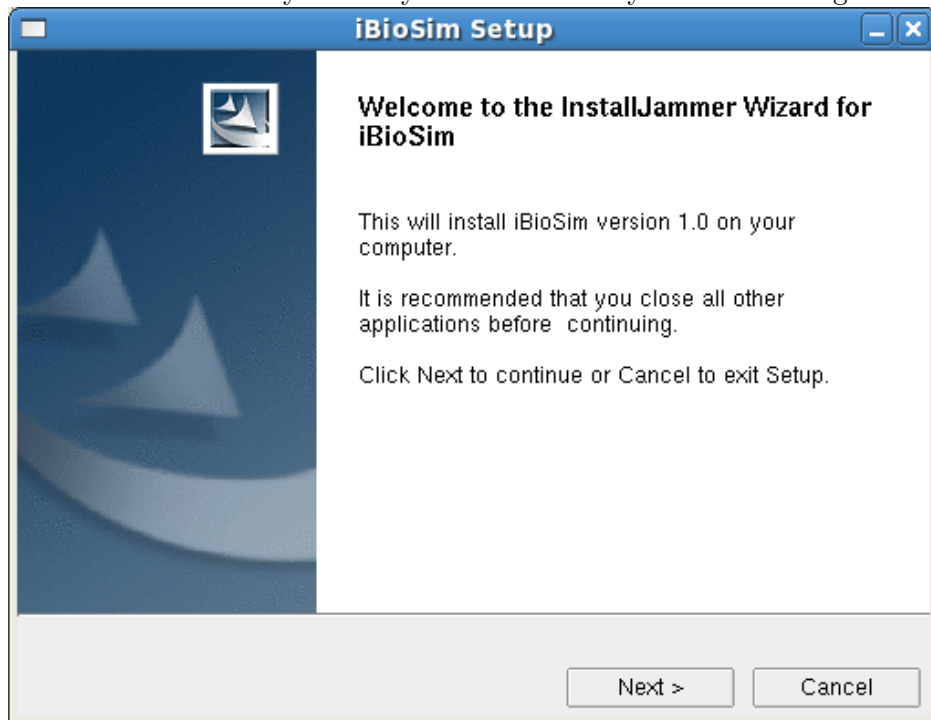


Next, it ask you for an installation location. Please make sure to select a path that does not have any spaces or special symbols as these cause problems with iBioSim.

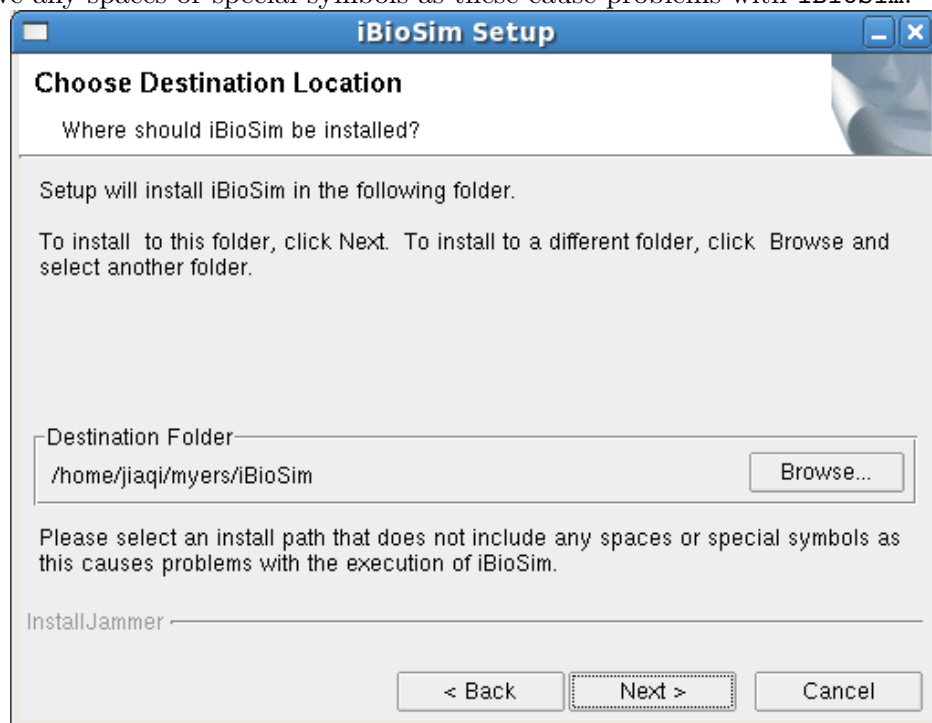

You are now ready to install. Press Next to continue.

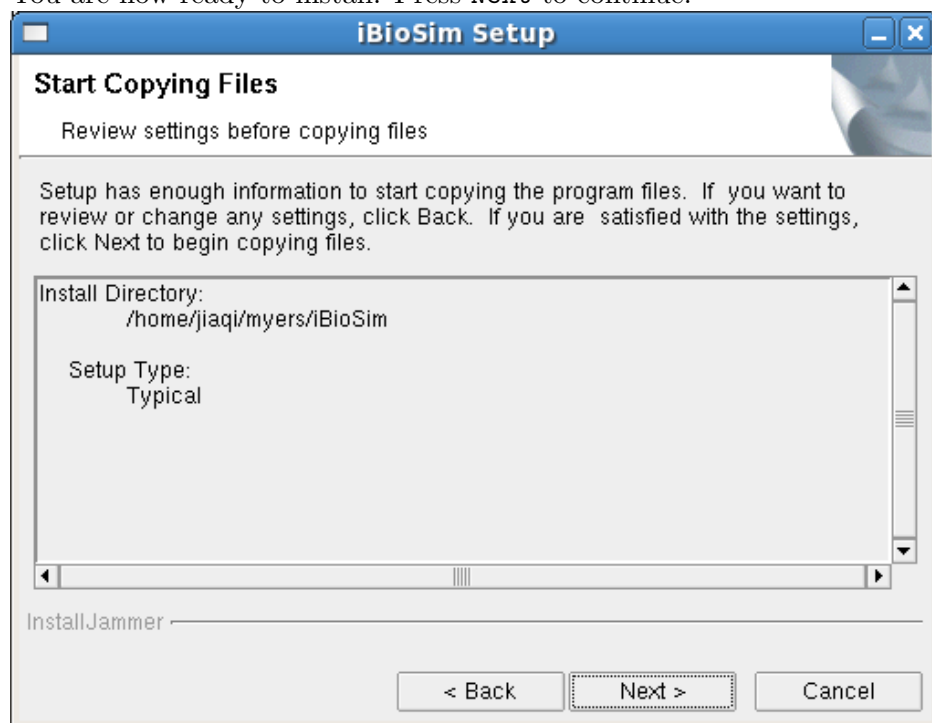

You are all done. Press **Finish**. If selected, iBioSim will launch immediately. Otherwise, you can start it using your desktop shortcut or from your start menu.

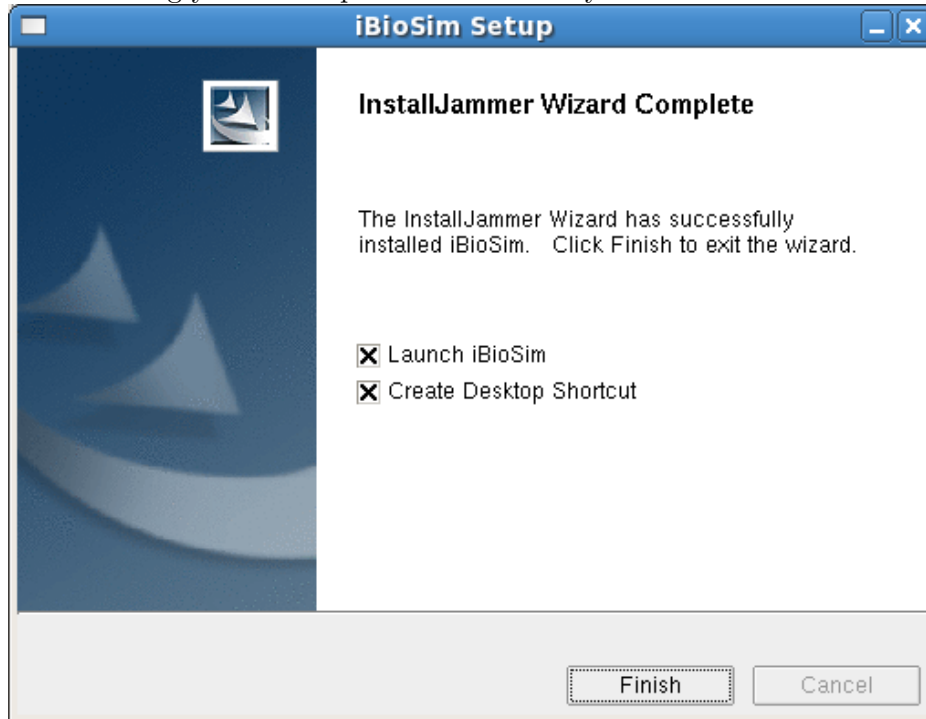

In your start menu, there should be an option to uninstall. If you select this, it will ask if you are sure then proceed to completely remove iBioSim from your system. It is highly recommended that you remove iBioSim using this uninstall procedure before installing a new version.

3 Installation on Linux


Since InstallJammer is also used for the Linux install, the installation instructions are essentially the same. First, download iBioSim-`<version>`-Linux-x86-Install. Open a terminal and browse to where this file was download. You must make this file executable:

`chmod u+x iBioSim-<version>-Linux-x86-Install`. You should then execute this file:


`./iBioSim-<version>-Linux-x86-Install`. This starts InstallJammer. It first asks you for your preferred installation language. Make your selection and press OK.




The next screen tells you that you what version you are installing. Press **Next** to continue.

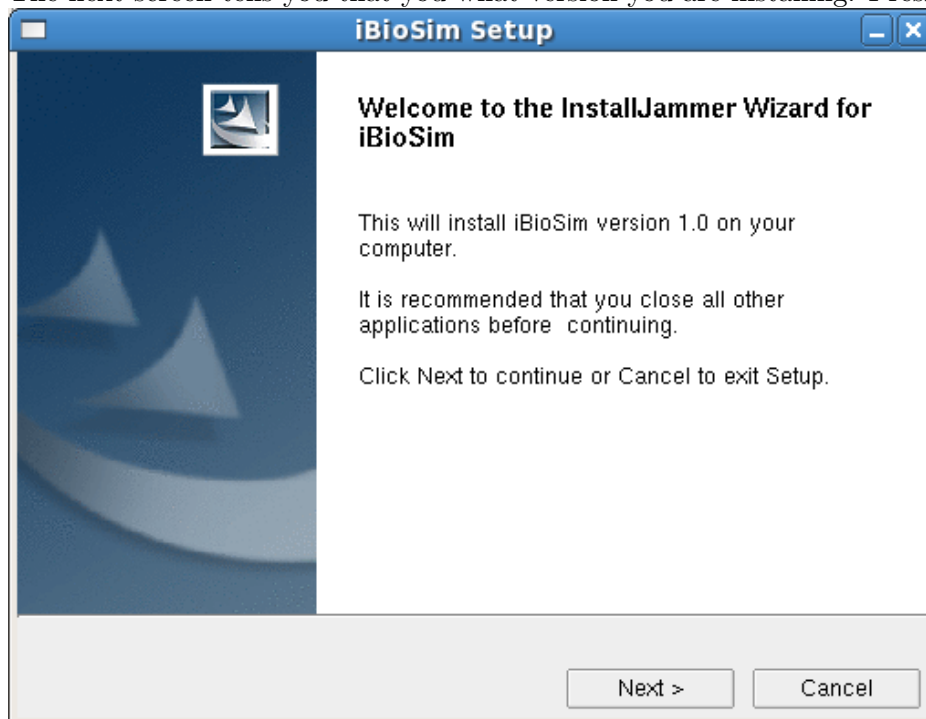


Next, it ask you for an installation location. Please make sure to select a path that does not have any spaces or special symbols as these cause problems with iBioSim.

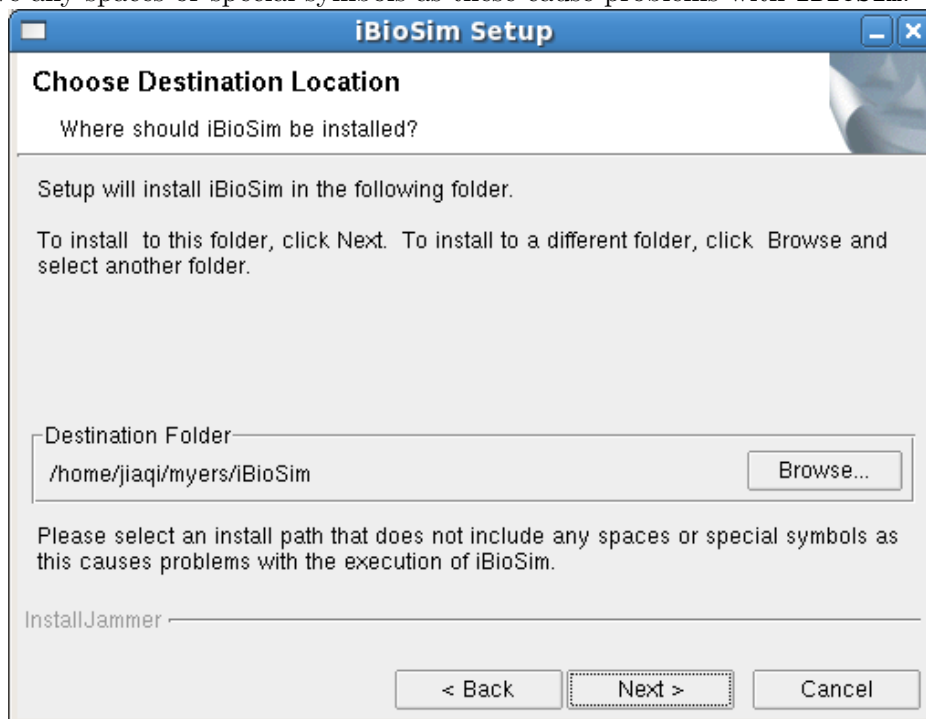

You are now ready to install. Press **Next** to continue.

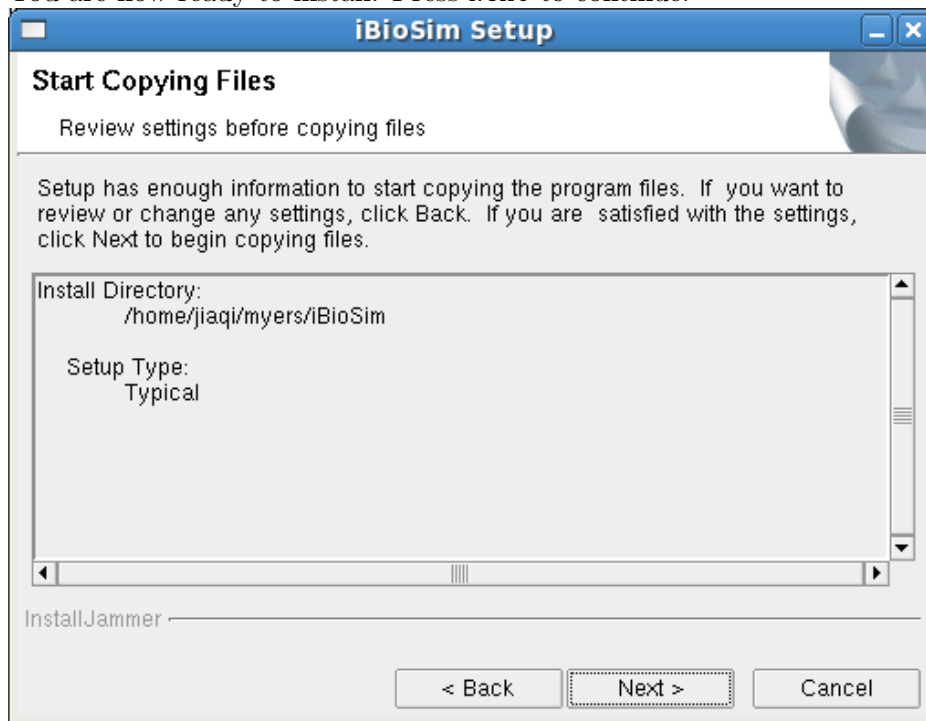

You are all done. Press **Finish**. If selected, iBioSim will launch immediately. Otherwise, you can start it using your desktop shortcut or from your start menu.

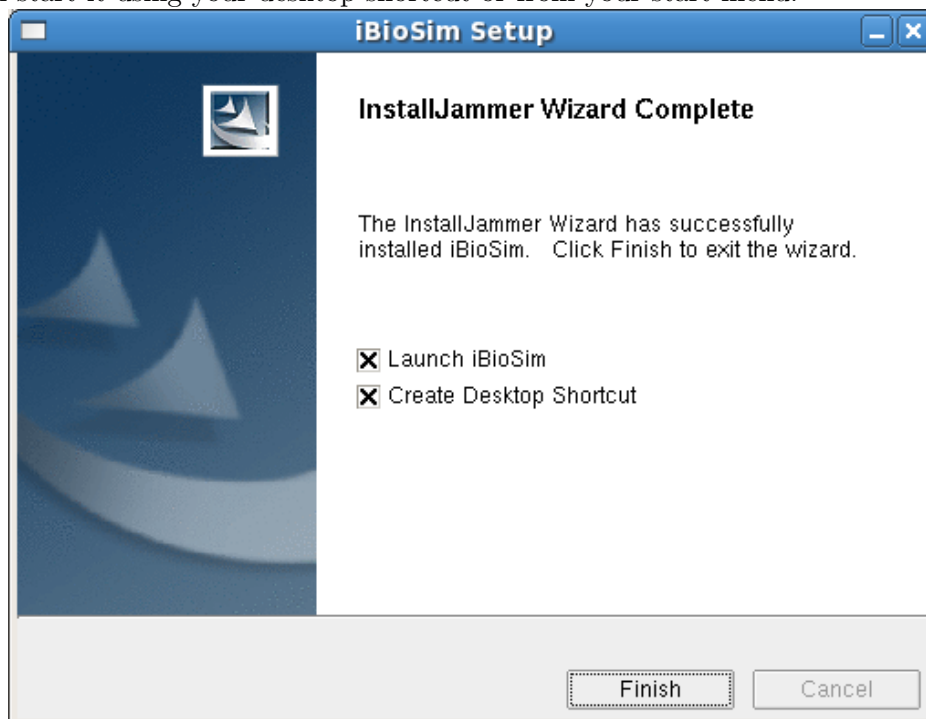

In your start menu, there should be an option to uninstall. If you select this, it will ask if you are sure then proceed to completely remove iBioSim from your system. It is highly recommended that you remove iBioSim using this uninstall procedure before installing a new version.

Note that since there are a large number of different linux installations, there can be problems

especially with older installations. If you are having problems though, first log out and back in. There are some environment variables that need to be set. If you are still having problems, then read on. Our linux installations are currently built for 32-bit and 64-bit versions of Fedora 15. While the Java parts of our program port smoothly, there are often issues with the dynamic libraries and the C/C++ programs **reb2sac** and **GeneNet**. If you are having a problem, the first step is to download and build libsbml Version 5.0. You will need to replace the libsbml.so* and libsbmlj.so files with your build of libsbml in either lib or lib64 depending on whether you are using a 32-bit or 64-bit machine. At this point, you should be able to start **iBioSim**. However, if you are still having problems, send an email to myers@ece.utah.edu including information about your OS and any error messages when you try to execute **iBioSim** from the command line.


You should also now test **reb2sac** and **GeneNet**. In your install directory, try the following:

- `cd bin`
- `export LD_LIBRARY_PATH=$BIOSIM/lib:$LD_LIBRARY_PATH`
Note that this is for bash shell, so you may need to change for your shell, if not bash.
- `./reb2sac`
- `./GeneNet`


You should see the usage lines for these programs. If you do not, you may need custom builds, contact myers@ece.utah.edu.

4 Installation on MacOS



To install on MacOS, you need to download `iBioSim.dmg`. You should open this file with `DiskImageMounter.app`. This should open up this disk image in `finder`. You should then double-click on the `iBioSim` package.

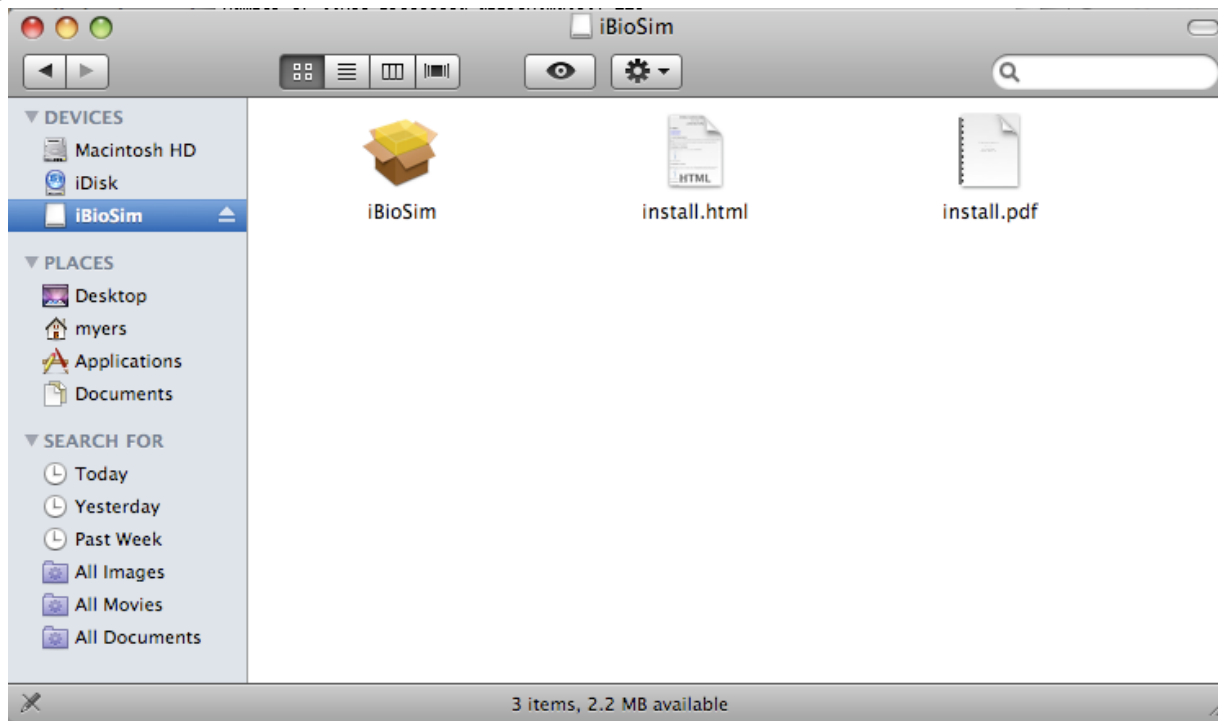


This starts the `iBioSim` installer. Press `Continue` to go on.

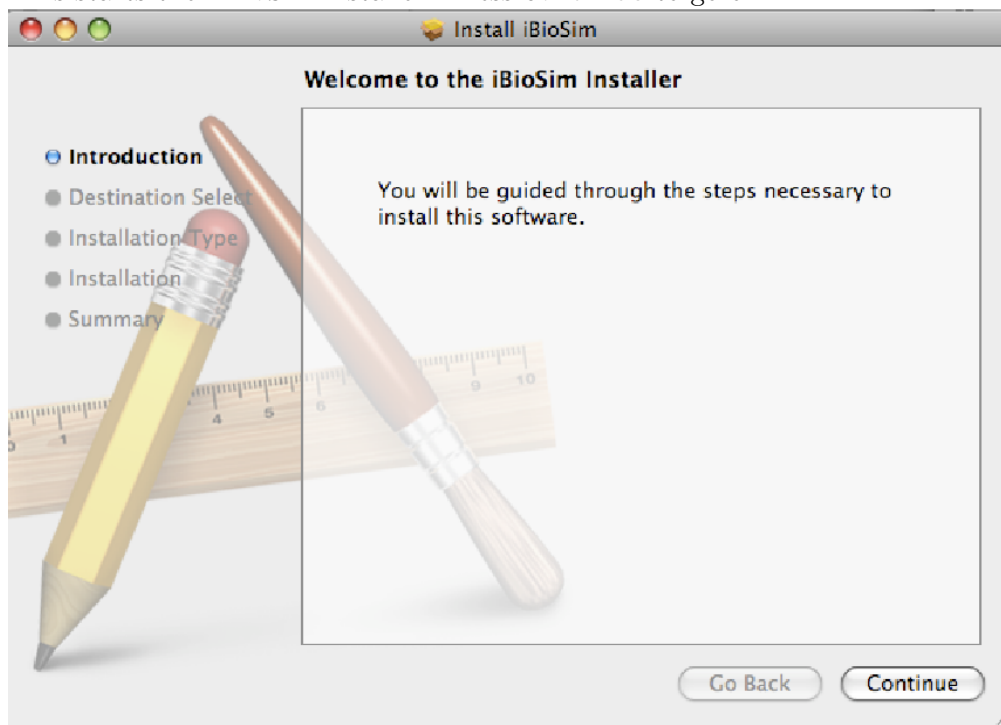

You should now select the drive you wish to install on. iBioSim will be installed in your Applications directory on this drive.

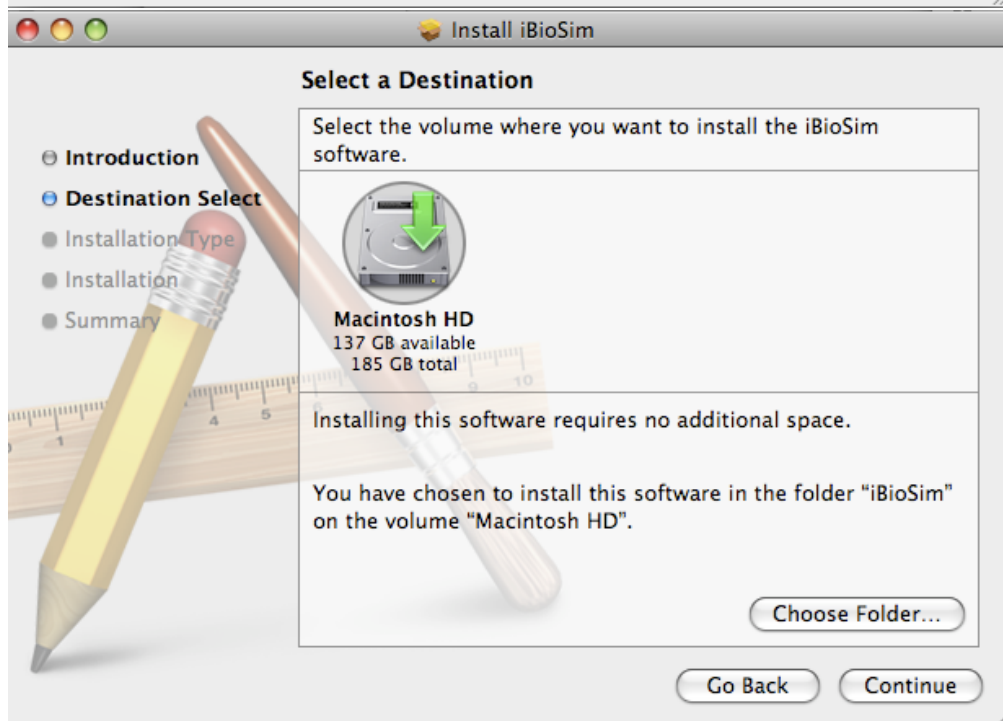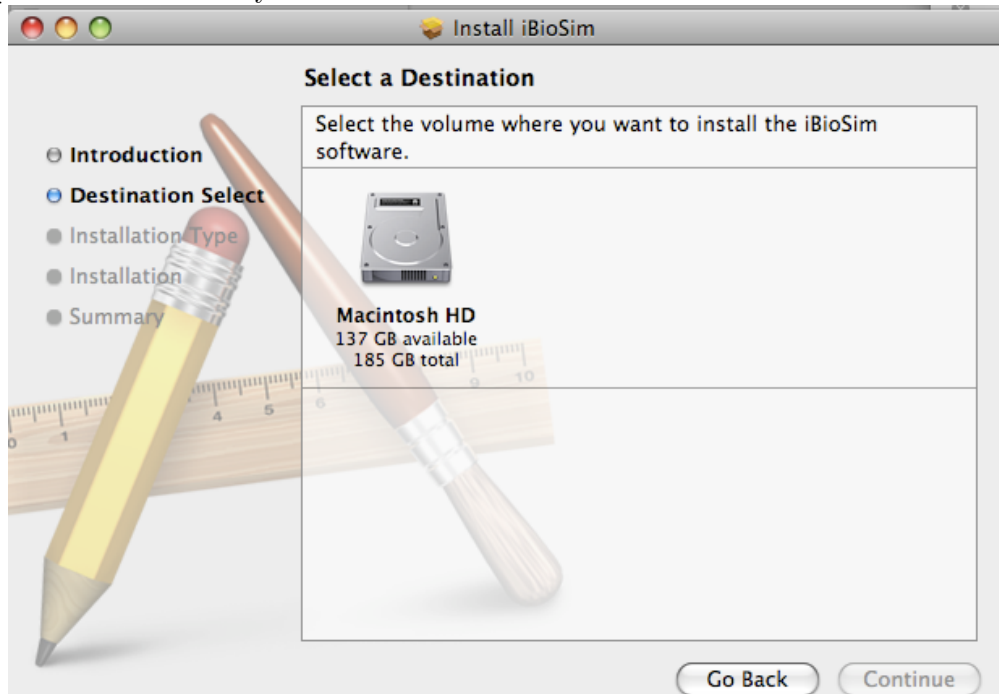

It will then tell you how much space it will take and allow you to confirm that you wish to install by pressing **Install**.

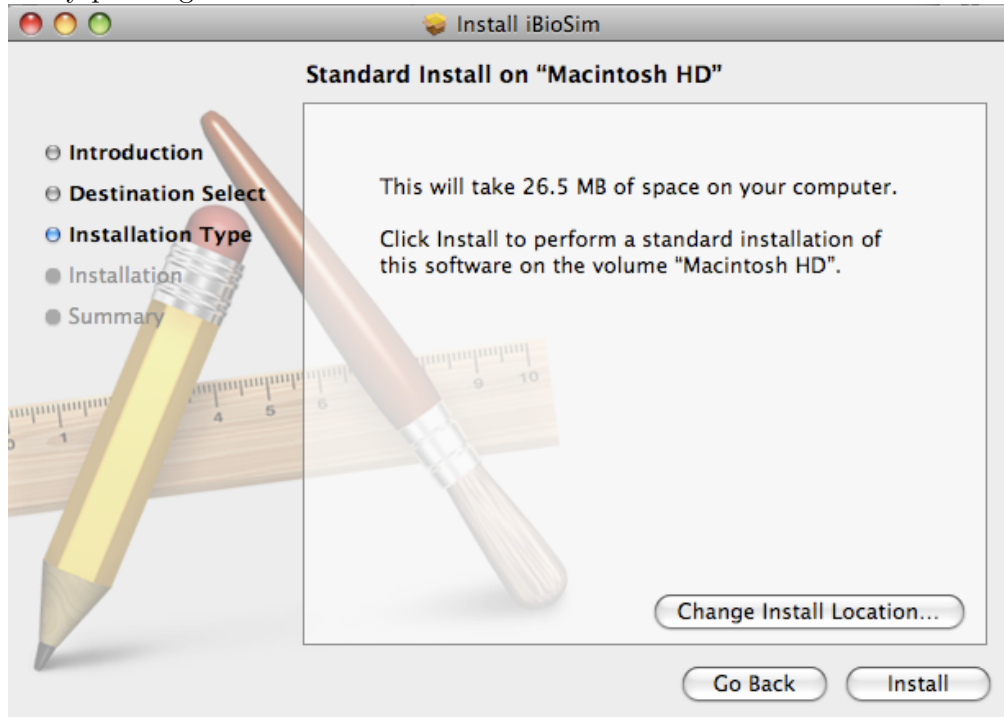

You should then get a message that it installed correct. At this point, you should be able to use Finder and locate iBioSim in your Applications directory. iBioSim is started by double-clicking on it.

